

BEFOUR

MEASURING NEW HEIGHTS

Display Features



- Automatically turns on when moved.
- Mounted on flexible neck for optimal viewing
- Selectable unit of measure
- Selectable decimal places



Made in the
USA

BEFOUR, INC


102 Progress Drive
Saukville, WI 53080
Phone: 800-367-7109
Fax: 262-284-5966
Email: sales@befour.com
www.befour.com



DIGITAL HEIGHT ROD: Model HTR-990 Digital Stadiometer Wall-Mount

HTR-990 Digital Stadiometer Wall Mount represents the ultimate digital stadiometer / height measurement system. The digital height rod provides unsurpassed 0.001 inch accuracy with choice of output in inches, centimeters or millimeters. The HTR-990 is a clear choice when you want the best in height measurement.

Slide mechanism rides on dual rail / quad bearing track for precise movement and is counterweighted just like an elevator for effortless one-finger travel. Digital stadiometer has best-in-class 62" selectable measurement span and features no button operation - simply move the slide to turn on. Display will automatically turn off after one minute with no activity. The display is mounted on a flexible neck for easy viewing regardless of the patient's height.

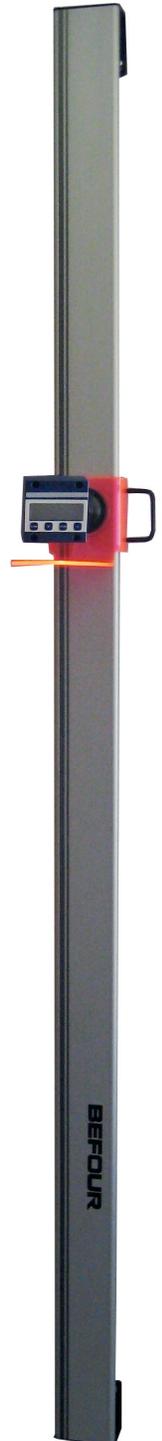

The measurement rail is constructed of attractive anodized aluminum for years of dependable use. The head-piece is constructed in safety-red acrylic with two convenient handles so that the digital stadiometer can be operated from either side. included are a 48" calibration level and all necessary mounting hardware.

FEATURES

- **Unsurpassed accuracy:** ± 0.001 inch
- **Choice of measurement unit:** inch, centimeters, millimeter
- **Best-in-Class measuring span:** Selectable 62" inch (i.e 24" to 86")
- **Precision slide mechanism:** Rides on dual/quad bearing system
- **Smooth, effortless travel:** Measurement slide has two convenient handles and is counterweighted. Remains stationary when not in use.

SPECIFICATIONS

Accuracy	± 0.001 in (0.003 cm)
Graduation (selectable)	
Inch	0.1, 0.01, 0.001, 1/16, 1/32, 1/64
Centimeter	0.1, 0.01, 0.001
Millimeter	0.1, 0.01
Measuring Span (selectable)	62 in (158 cm)
Display	0.5 in LCD
Power	(1) CR123 Lithium Battery
Rail Dimensions (HxWxD).....	76" x 4.5"x 12" (193 x 11 x 31 cm)
Country of Origin	Made in the USA
Warranty	3 years