

Professional Healthcare Scales

BEFOUR

1-800-367-7109

www.befour.com


CX7 CONSOLE FEATURES:



- Displays weight, height, and BMI on a single screen
- EMR connectivity via serial output or optional wireless module
- ECO Weigh battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)
- 3-segment battery life indicator
- Large 1" LCD digits
- No AC outlet or recharging needed



Befour products are made with pride in the USA!

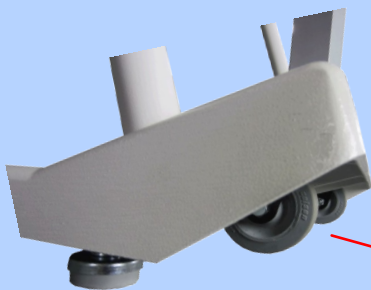

To move to a new location, simply tilt the scale back on its large 3" wheels and roll away. The scale fits through all standard doorways.

Model **MX310**

Tilt & Roll Handrail Scale with BMI

The finest handrail scale on the market, it is ideal for most patient weighing requirements. With a 1000 pound capacity and ergonomic handrails, the patient's weight is determined accurately, quickly, and with maximum safety. The "tilt & roll" design, large wheels, and battery operation makes moving the scale effortless. Operation is easy and intuitive, even for the first time user.

The standard, integrated handrail can be grasped by patients while still providing an accurate measurement. The scale features Befour's proprietary MotionLock Technology which automatically calculates the weight of the most unsteady patients. It's alkaline battery power provides years of use and allows the scale to

FEATURES/SPECIFICATIONS

Wrap-around heavy-duty handrail: Added security for patients of all sizes. Unique handrail design does not interfere with foot positioning.

Fully assembled: Arrives ready to use, just put in the leveling feet and batteries.

"Tilt & Roll" Design: Always ready to move to another location.

MotionLock Technology: Accurate weight even with an unsteady user.

LB/KG Lockout: Quickly change setup to disable LB/KG button so scale will only weigh in selected unit of measure.

New CX7 console: Displays only the information you need - weight only or weight, height, and BMI.



Ergonomic, integrated handrails allow patient to grasp anywhere for support. Patient's weight is accurately calculated even when using the handrails.



Scale Height	42"
Platform Dimensions	22" x 25" x 2"
System Weight	46 lbs.
Weight Capacity	1000 lbs. / 450 kg
Resolution	0.1 lb. / (0.05 kg.)
Accuracy	± 0.1 lb. / (0.05 kg.)
Power Source	4 D-cell alkaline batteries
Connectivity.....	USB Serial or optional wireless
Warranty	3 Years

BEFOUR