




BARIATRIC SCALE

High-capacity, portable Pro-BMI scales are designed with precision and heavy weighing in mind.

Proudly Made in the USA

Model PS-7700 Pro-BMI Scale

Constructed with oversized platform and triple-reinforced undercarriages...

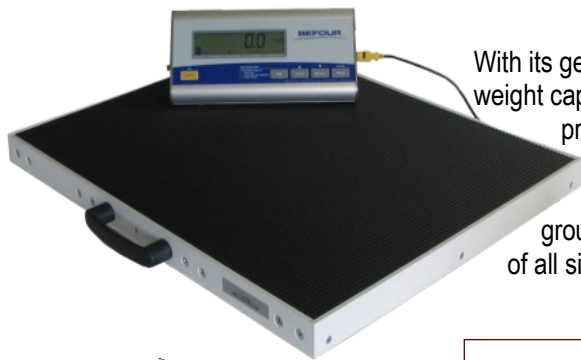

With its generous 24" x 20" platform and a 900 lb. weight capacity, the PS-7700 is ideal for bariatric professionals, home health care, weight-lifting competitions, and athletic weigh-ins. The low-profile platform sits just 2" off the ground, making it easy to step on for individuals of all sizes.



MODEL PS-7700

Platform Dimensions.....	20" x 24" x 2"
Weight Capacity.....	900 lbs. or 405 kg.
Accuracy / Resolution.....	± 0.1 lb. (± 0.1 kg.)
Power Source.....	6 C batteries
Readout.....	1" LCD digits
System Weight.....	25 lbs.
Warranty.....	3 Years

FEATURES/SPECIFICATIONS:

Body Mass Index (BMI) Calculator: Calculates body mass index in seconds in two easy-to-use modes. Standard Mode, used in most situations.

Assistant Mode, allows a fitness trainer or health professional to enter a person's height before weighing.

Unique MotionLock Technology: Incredibly accurate measurements even for unsteady users.

High Capacity Tare: Resets scale to zero, even with weight on platform (i.e. mother-child function).

Push Button Calibration: Set calibration via keypad entry without the need of tools.

Programmable Set Up: Tailor the console for your specific needs including auto-off time, audible beep, lb. or kg. lockout, and more.

Lightweight and Easy to Carry: Built-in handle and recessed storage cavity securely holds the readout during storage and extended transport.

Weight Measurements: In pounds or kilograms.



PERFECT FOR WEIGHING
THE ELDERLY
&
HEAVY SET!

LCD DISPLAY FEATURES:



- Easy to use Body Mass Index (BMI) calculation with programmable clothing allowance
- Pinpoint accuracy of ± 0.1 lb.
- The console is powered by six (6) standard C-cell alkaline batteries (included)
- Unique MotionLock Technology has incredibly accurate measurement even for the most unsteady user
- Large 1" LCD digits
- Auto-off time from 12 seconds to 25 minutes, audible tone when weight locks in, lb/kg lockout mode

