



Proudly Made in the USA

HEALTH & FITNESS SCALE

Quick, reliable, and extremely accurate in high-volume environments...

Model PS-8070 Tilt & Roll Handrail Scale

Stand-on scales are ideal for locker rooms, fitness clubs, medical facilities, and personal use.

Rugged construction, portability...

The PS-8070 provides a truly versatile scale that is ideal for most weighing needs. It is equally at home at the health club, in the doctor's office, or on the shop floor - weighing people or packages. Convenient battery power allows user to place scale anywhere, no outlet or re-charging needed! Sturdy, integrated handrail is a standard feature that can be grasped by the user while still providing an accurate measurement. Transport is easy ... simply tilt the scale back on its heavy-duty wheels and roll to the next location.

- Quick 2 second weight measurement
- Extreme accuracy of ± 0.1 lb./0.1 kg
- Heavy-duty 1 1/4" diameter handrail provides added security
- Large 3" wheels roll effortlessly on all flooring surfaces
- Advanced design results in a scale of only 25 lbs.



MOTIONLOCK

FEATURES/SPECIFICATIONS:

Programmable Settings: Includes auto-off time, "always-on" option, audible tone, lb. or kg. lockout, and more.

Dual Tare Operation: Automatic and keypad tare features help determine net weight.

Auto-zero, weight recall and reweigh.


MotionLock Technology: Accurate weight displays even with unsteady users.

Push Button Calibration: Set calibration via keypad entry with no tools needed.

3 Year Warranty: An industry best.



Weigh People

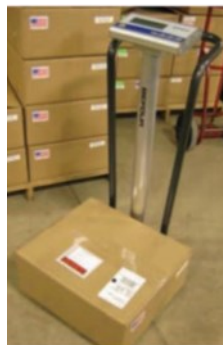

Weigh Packages

LCD DISPLAY FEATURES:



- Alkaline batteries provide over 80,000 weighings
- "Always-On" option
- Weight reads in pounds or kilograms
- Programmable pound/kilogram lockout
- Programmable auto-off: 12 seconds to 25 minutes
- Auto zero
- Weight recall and reweigh
- Audible sound indicates weight is "locked-in"
- Durable anodized aluminum construction
- Large 1" easy-to read digits
- 3-segment battery life indicator

MODEL PS-8070

Scale Height.....	42"
Platform Dimensions.....	16" x 18" x 1 3/4"
Weight Capacity.....	500 lbs. or 225 kg.
Accuracy / Resolution.....	± 0.1 lb. (± 0.1 kg.)
Power Source.....	6 C batteries
Readout.....	1" LCD digits
System Weight.....	25 lbs.
Warranty.....	3 Years



Tilt & Roll

Befour Inc.

102 Progress Drive
Saukville, WI 53080

Phone

800-367-7109

Web

www.befour.com

Fax

262-284-5966

E-mail

sales@befour.com