


Befour's HTR101 digital height rod features fast & accurate performance.

With its user-selectable 51" measurement span, the versatile HTR101 can measure the height of a small child to that of an NBA player. The smooth glide headpiece features one-button operation, bright LED digits that can be easily seen, even at maximum height, and a unique auto-hold function which retains the measurement display even as it is raised out of the way.

The HTR101 has an accuracy of ± 0.1 in (0.1 cm), easily converts between three display formats (inches, feet & inches or centimeters) and incorporates an automatic turn-off. The efficient battery system provides 10,000 height measurements on a single set of (2) "AAA" alkaline batteries (included).

With its seamless design and no visible wires, the HTR-101 is ideal for hospitals, clinics, educational institutions and fitness centers.

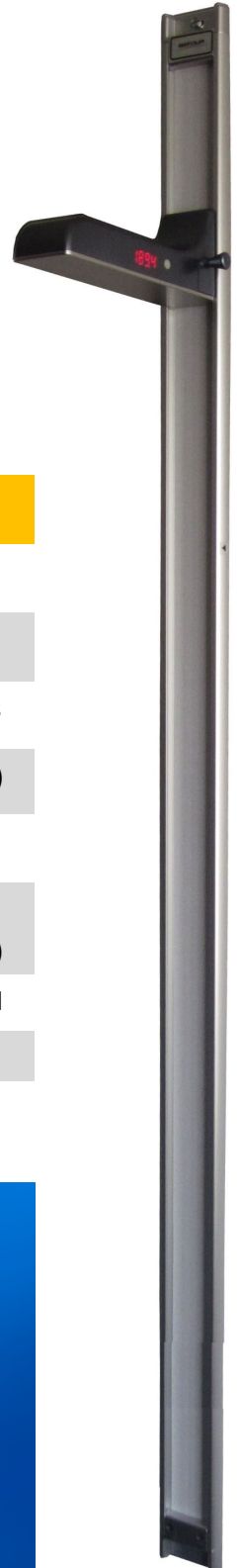

- Fast & accurate measurement
- View height in selectable formats
- Measures up to 91"
(7 ft 7 in / 230 cm)
- Effortless sliding headpiece
- Flexible measurement range



HEIGHT ROD SPECIFICATIONS

ACCURACY ± 0.1 in (0.3 cm)

MEASURING SPAN 51 in (130 cm)

DISPLAY Bright LED Digits

POWER (2) "AAA" batteries (included)

WEIGHT 2.5 lbs

RAIL DIMENSIONS (HxWxD) 57.5 x 2.2 x 0.5 in
(146 x 5.6 x 1.3 cm)

HEADPIECE Extends 8 inches from rail

CONSTRUCTION Anodized aluminum rail

WARRANTY 3 years

DISPLAY FEATURES

PICK YOUR DISPLAY FORMAT:

Measurements are displayed in BRIGHT LED digits in
Inches, Feet & Inches, or Centimeters



Inches



Feet & Inches



Centimeters