

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Portable Bariatric Scale with BMI

MODEL MX170

The ultra-rugged, oversized, MX170 scale - with 1,000 pound capacity - is ideal for bariatric clinics.

This scale has no detectible flex underfoot and the rock-solid platform incorporates triple-walled extruded aluminum framework and full-perimeter tensile reinforcement. The **extra-large platform** provides stability, while the built-in handle, battery power and display console storage cavity provide maximum portability. Operation is simple and intuitive, even for the first time user.

The scale offers many **advanced features which come standard**, including Befour's proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an **'industry leading' 3-year warranty** and are **Made in the USA**.

FEATURES & BENEFITS

Pin point accuracy of ± 0.1 lb.

High Capacity: 1,000 pounds

Fully Assembled: Arrives ready for use - simply insert the provided leveling feet and batteries


Portability: Built-in handle and recessed storage cavity for display console, create an easy-to-transport weighing system

Built to Last: Rugged, non-corroding all-aluminum design weighing only 25 pounds

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required.

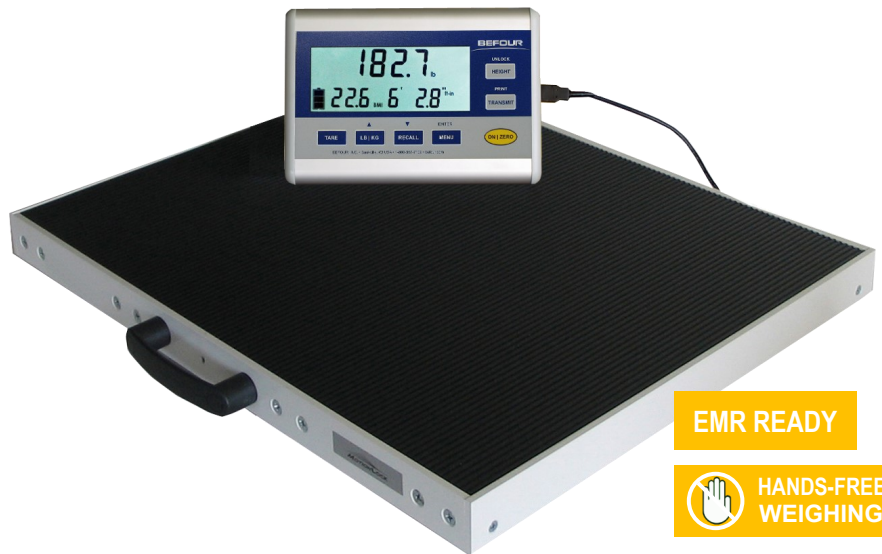

EMR READY



HANDS-FREE
WEIGHING



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or optional wireless module

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

Portable Bariatric Scale with BMI

MODEL MX170

SCALE SPECIFICATIONS

CONSOLE MOUNTING Wall or tabletop

PLATFORM DIMENSIONS 24" W x 20" D x 2" H

SYSTEM WEIGHT 25 lbs

WEIGHT CAPACITY 1000 lbs (450 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available (P/N: 03049-08)

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or optional wireless

DISPLAY READOUT 1" LCD digits

WARRANTY 3 years

SHIPPING SPECIFICATIONS


PACKAGED WEIGHT 31 lbs

PACKAGE DIMENSIONS Box (30" x 25" x 4")

SHIPPING METHOD UPS Package

DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

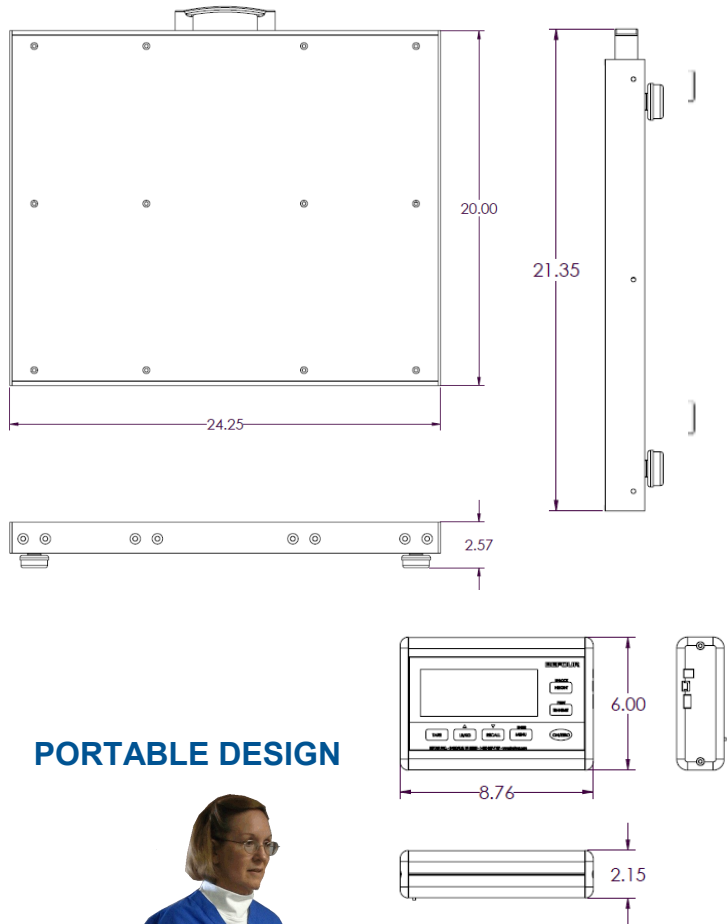

PORTABLE DESIGN



Versatile Console:

Mounting brackets (attached) allow positioning on wall or angled on tabletop



Display console can be securely stored underneath scale platform and battery operation eliminates need for wall outlets, making this scale extremely portable.

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com