

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Fold & Roll Wheelchair Scale with Handrail

MODEL MX480 & MX480D

This oversized, heavy-duty scale can be used as a wheelchair, handrail, or chair scale. Single ramp can be mounted on either side, while the dual ramp version allows for pass-through weighing. The ramps also fold up for use as a stand-on scale.

The **“fold & roll” design and large wheels** makes transportation of the scale effortless. Operation is simple and intuitive, even for the first time user.

The scale offers many **advanced features which come standard**, including Befour’s proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an ‘industry leading’ 3-year warranty and are **Made in the USA**.

FEATURES & BENEFITS

Generous platform size: 32” x 34” platform accommodates most wheelchairs

Pivoting Console: Pivots 320° for optimal viewing

Folding Design: Platform folds up to take up less space

Integrated “Active” Handrail: Allows patient to grasp handrail for support while weight is accurately calculated


Fully Assembled: Arrives ready for use; simply insert the provided leveling feet and batteries

“Fold & Roll” Design: Safely weigh patient when scale is stationary, but ready to move at any time

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required

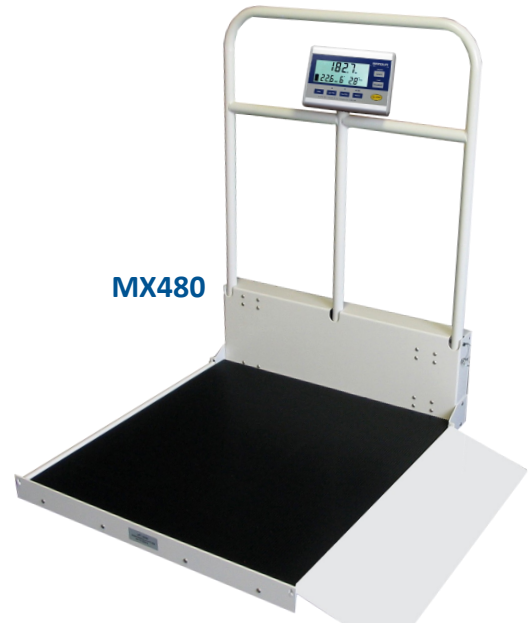

MX480



MX480D

EMR READY



HANDS-FREE
WEIGHING



DISPLAY CONSOLE

Display is used across all Befour’s medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output of other options

Programmable User Settings can be set to auto-off, “Hands-Free” Weighing, etc.

‘ECO Weigh’ battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1” LCD digits

No AC outlet or recharging needed

Fold & Roll Wheelchair Scale with Handrail

MODEL MX480 & MX480D


SCALE SPECIFICATIONS

SCALE HEIGHT	45"
PLATFORM DIMENSIONS	32" W x 34" D x 2" H
SYSTEM WEIGHT	87 lbs
WEIGHT CAPACITY	1000 lbs (450 kg)
RESOLUTION	0.1 lb (0.05 kg)
ACCURACY	± 0.1 lb (0.05 kg)
POWER SOURCE	4 D-cell alkaline batteries (incl) <small>*optional AC adapter available (P/N: 03049-08)</small>
BATTERY LIFE	100,000 weighings
CONNECTIVITY	USB serial or other options
DISPLAY READOUT	1" LCD digits
WARRANTY	3 years
PACKAGED WEIGHT	103 lbs
PACKAGE DIMENSIONS	Box (40" x 35" x 14")
SHIPPING METHOD	Freight


DIMENSIONS

Allow 1" perimeter space - scale should not touch wall


MX480 = Single Ramp

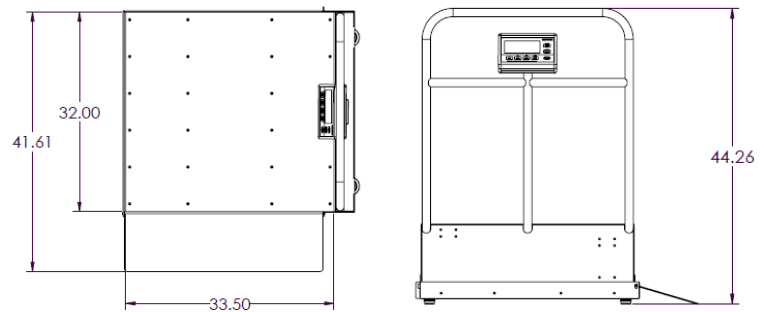

MX480D = Dual Ramp

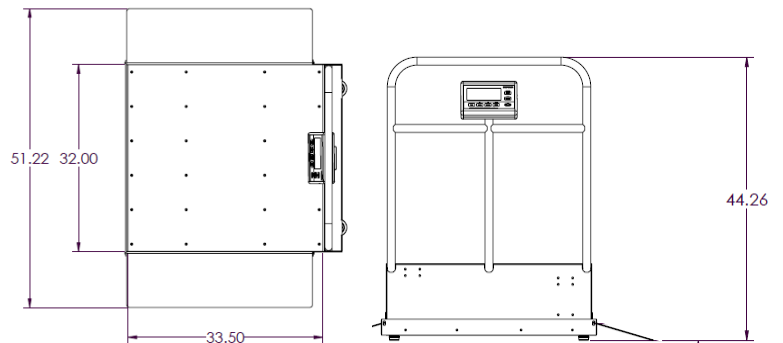

SHIPPING SPECIFICATIONS

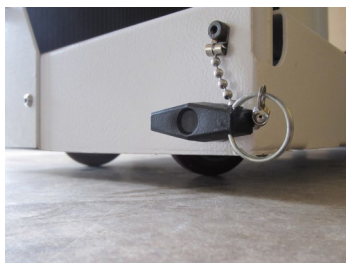

Folded scale rolls effortlessly to the next location on four hidden wheels, two - pivot.

FOLD & ROLL TRANSPORT



Footprint of folded scale is 32" x 12", is self-supporting for safe transport and storage, and will fit through all standard doorways.

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com