

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Tilt & Roll Handrail Scale with BMI & Height Rod

MODEL MX810

The MX810 is the premier mobile measuring station on the market. With its integrated digital height rod, this scale quickly determines weight, height & BMI. It features a 1,000 pound capacity and an ergonomic handrail that can be grasped during weighing.

The “tilt & roll” design makes transportation of the scale effortless. Operation is simple and intuitive, even for the first time user.

The scale offers many **advanced features which come standard**, including Befour’s proprietary MotionLock technology and ECO Weigh battery system, for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an ‘industry leading’ 3-year warranty and **USA-Made**.

FEATURES & BENEFITS

Integrated “Active” Handrail: Allows patient to grasp handrail for support while weight is accurately calculated

Ergonomic Wrap-Around Handrail: Unique, heavy-duty, design does not interfere with foot positioning and adds security and stability for patients of all sizes

Fully Assembled: Arrives ready for use - simply insert the provided leveling feet and batteries

“Tilt & Roll” Design: Safely weigh patient when scale is stationary, but ready to move at any time

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more

Digital Height Measurement: Simply press HEIGHT & raise the foldable headpiece. Height can be viewed in selectable units & displays automatically



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required



DISPLAY CONSOLE

Display is used across all Befour’s medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, “Hands-Free” Weighing, etc.

‘ECO Weigh’ battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1” LCD digits

No AC outlet or recharging needed

Tilt & Roll Handrail Scale with BMI & Height Rod

MODEL MX810

SCALE SPECIFICATIONS


SCALE HEIGHT	42"
PLATFORM DIMENSIONS	22" W x 25" D x 2" H
SYSTEM WEIGHT	48 lbs
WEIGHT CAPACITY	1000 lbs (450 kg)
RESOLUTION	0.1 lb (0.05 kg)
ACCURACY	± 0.1 lb (0.05 kg)
POWER SOURCE	4 D-cell alkaline batteries (incl) <small>*optional AC adapter available (P/N: 03049-08)</small>
BATTERY LIFE	100,000 weighings
CONNECTIVITY	USB serial or other options
DISPLAY READOUT	1" LCD digits
HEIGHT ROD RANGE	38-88" / (97-225cm)
HEIGHT ROD ACCURACY	± 0.1 in (0.2 cm)
HEIGHT DISPLAY FORMATS	In, Ft & In, Cm
WARRANTY	3 years

SHIPPING SPECIFICATIONS

PACKAGED WEIGHT	70 lbs
PACKAGE DIMENSIONS	Box (27" x 24" x 43")
SHIPPING METHOD	Freight

DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

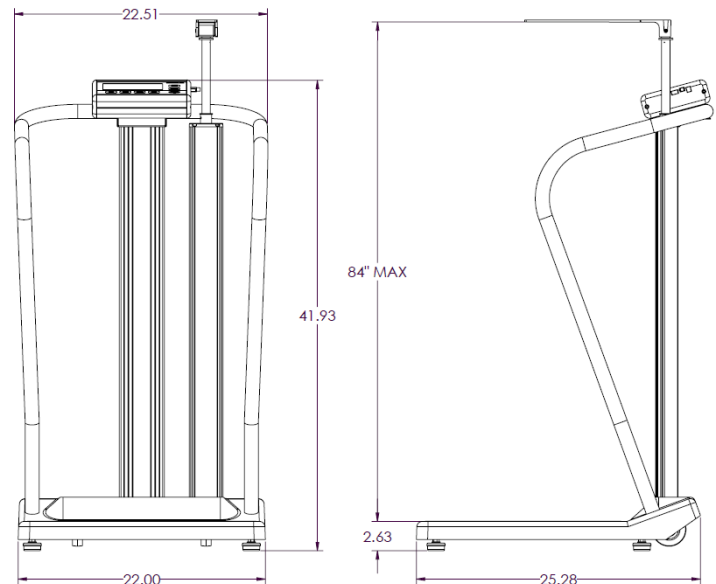

To move to a new location, simply tilt the scale back on its large 3" wheels and roll away. Scale fits through all standard doorways.



Height rod is both digital and integrated into scale for automatic BMI calculation. Unique design has extended measuring range with no exposed cables and rotates 360° for easier positioning of headpiece.

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com