

BEFOUR

PROFESSIONAL HEALTHCARE SCALES


Wall Mount 'Fold Away' Wheelchair Scale

MODEL MX440D

The MX440D Wall Mount Scale is a space-saving unit that adheres to the wall and tucks away neatly when not in use to maximize floor space. Integrated struts allow the extra-large weighing platform to lower with ease and raise back to folded position with soft close technology. The integrated dual ramps simply fold up and down, making this scale a versatile option for either wheelchair or stand-on weighing, while accommodating patients with walkers or canes.

This scale is a 'fold away' design and offers many **advanced features which come standard**, including Befour's proprietary MotionLock technology and ECO Weigh battery system allowing 100,000+ weight readings on a single set of batteries.

All Befour scales are **built to last**, constructed of heavy-duty metal, backed by an **'industry leading' 3-year warranty** and are **Made in the USA**.

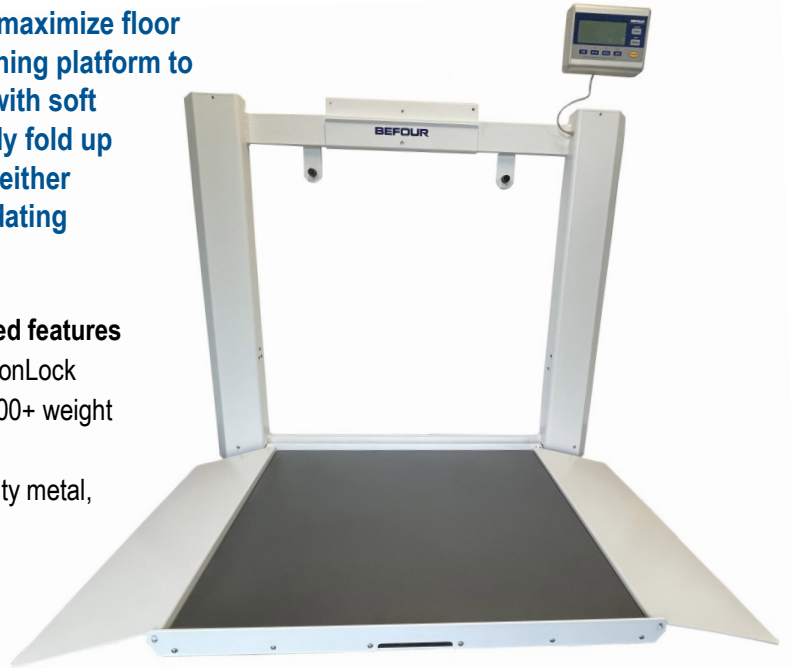

FEATURES & BENEFITS:

Generous platform size: 32" deep x 32" wide platform accommodates wheelchairs and gurneys

Space Saving: Tucks neatly against wall when not in use

Lift Assist: Actuators make raising and lowering the platform a breeze

Simple Assembly: Minimal set-up required

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required

EMR READY



HANDS-FREE
WEIGHING



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or optional wireless module

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed


Wall Mount 'Fold Away' Wheelchair Scale

MODEL MX440D

SCALE SPECIFICATIONS

SCALE HEIGHT	42"
PLATFORM DIMENSIONS	32"W x 32"D x 2"H
SYSTEM WEIGHT	90 lbs
WEIGHT CAPACITY	1000 lbs (450 kg)
RESOLUTION	0.1 lb (0.05 kg)
ACCURACY	± 0.1 lb (0.05 kg)
POWER SOURCE	4 D-cell alkaline batteries (incl) *optional AC adapter available (P/N: 03049-08)
BATTERY LIFE	100,000 weighings
CONNECTIVITY	USB serial or optional wireless
DISPLAY READOUT	1" LCD digits
WARRANTY	3 years
SHIPPING SPECIFICATIONS	
PACKAGED WEIGHT	~145 lbs
PACKAGE DIMENSIONS	Box approx. (48" x 52" x 7")
SHIPPING METHOD	Freight

DIMENSIONS

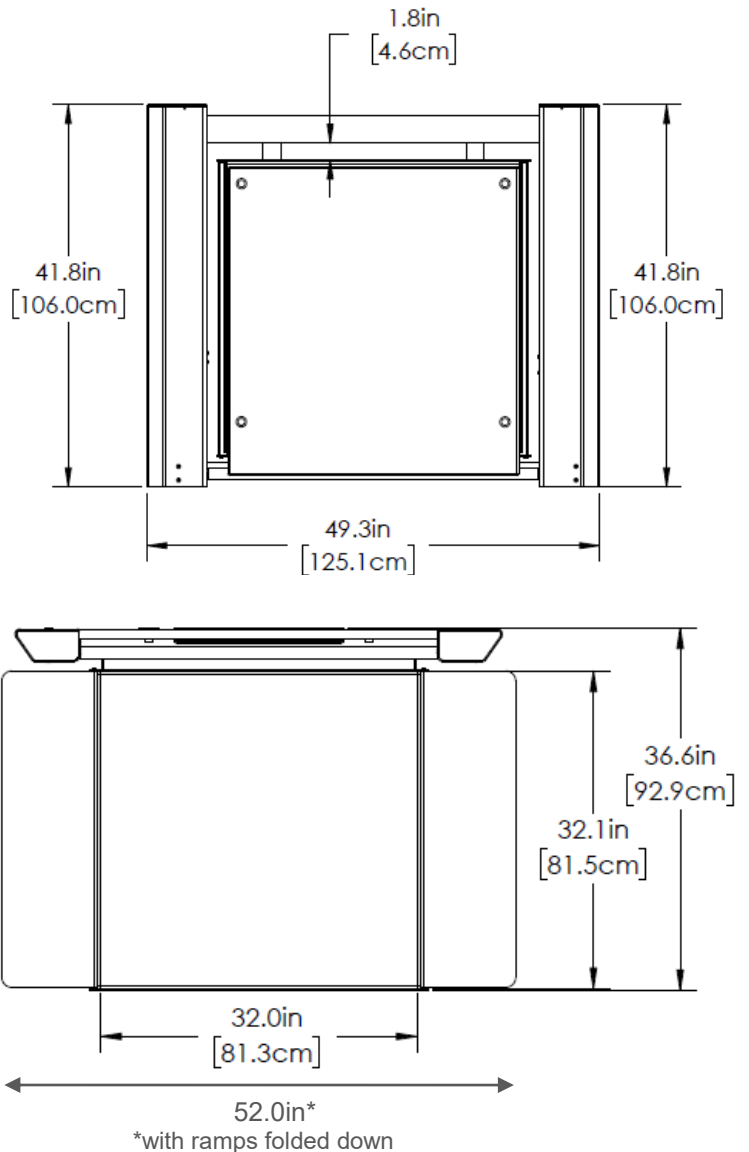

**Folded Position
Against Wall**

Ramps fold down when in use



BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com